



VitalJacket® SDK


VitalJacket SDK v1.0.07 – QRS detector

LEGAL NOTICE AND DISCLAIMER


ATTENTION: Although **VitalJacket** is a certified medical device, its developer version is **NOT** certified for diagnosis usage. It is intended for R&D and development purposes **only**. Users of VJ SDK can submit their final developments to medical certification. All contents of our product are compliant with the European Medical Device directive 93/42/EEC but, being a developer's version, it's not certified.

Index:

QRS detector	4
Class: Detector QRS	5
Class: QRS info	7
Sample Code	8
References	10
Control versions	11

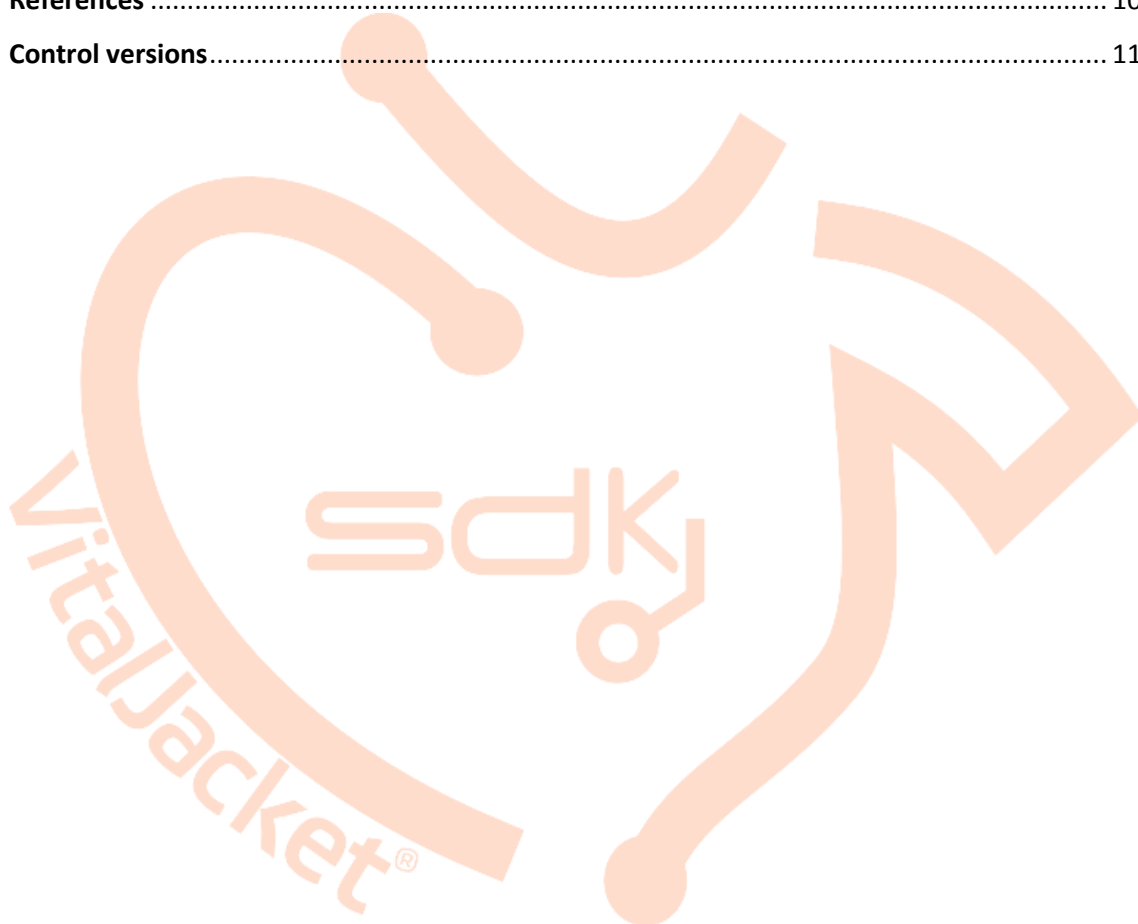

QRS detector

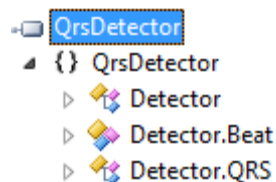

public class **Detector**
Member of [QrsDetector](#)

Summary:

QRS detector. The QRS detector is based on the algorithm of Pan and Tompkins [1] and was used MIT-BIH database for validate results [2].

public class **QRS**
Member of [QrsDetector.Detector](#)

Summary:

Class: QRS.

Class: Detector QRS

- Detector(int, float, float, byte)
- GetPeak()
- GetPulse(float)
- GetQRS()
- GetRRR(int, int)
- QRSDet(int, short, int)

public **Detector**([int](#) sampleFrequency, [float](#) thPeakAmplitude, [float](#) th, [byte](#) windowSize)
Member of [QrsDetector.Detector](#)

Summary:
Construtor.

Parameters:
sampleFrequency: Sample frequency (Hz)
thPeakAmplitude: Threshold of peak amplitude
th: Threshold for detect noise.
windowSize: Number of QRS for calculate mean of Heart Rate

public [short](#) **QRSDet**([int](#) datum, [short](#) init, [int](#) sampleCount)
Member of [QrsDetector.Detector](#)

Summary:
Detect QRS.

Parameters:
datum: ecg sample
init: '1' for init QRS detector and '0' for detect QRS
ind: number of ecg samples

public [QrsDetector.Detector.QRS](#) **GetQRS**()
Member of [QrsDetector.Detector](#)

Summary:
Get last QRS detected.

Returns:
QRS detected

public [QrsDetector.Detector.Beat](#) GetPeak()
Member of [QrsDetector.Detector](#)

Summary:

Get last peak detected.

public static [float](#) GetPulse([float](#) rr)
Member of [QrsDetector.Detector](#)

Summary:

Calculate bpm instantaneous (bpm).

Parameters:

rr: R-R (ms)

Returns:

value bpm

public static [float](#) GetRR([int](#) rr, [int](#) SampleFrequency)
Member of [QrsDetector.Detector](#)

Summary:

Calculate R-R (ms).

Parameters:


rr: R-R in samples

SampleFrequency: Sample frequency (Hz)

Returns:

R-R (ms)

Class: QRS info

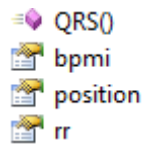

public [int](#) bpm { set; get; }
Member of [QrsDetector.Detector.QRS](#)

Summary:
Bpm (bpm).

public [int](#) position { set; get; }
Member of [QrsDetector.Detector.QRS](#)

Summary:
Peak position (samples).

public [int](#) rr { set; get; }
Member of [QrsDetector.Detector.QRS](#)

Summary:
R-R (ms).

Sample Code

1. Add Reference QrsDetector.dll to your project;
2. See sample code:

```
private int sampleCount = 0, windowCounter = 0;
private int sampleFrequency = 500;
private Queue<int> pulseQueue = new Queue<int>();
private int oldPeak = 0;
private float pulseValue = 0;
private Detector detector;

private void InitQrsDetector()
{
    detector = new Detector(500, 7.0f, 0.3125f, 8);
}

/// <summary>
/// Detect QRS.
/// </summary>
/// <param name="dataByte">ecg sample</param>
private void QrsDetector(short dataByte)
{
    try
    {
        sampleCount++;

        if (windowCounter > sampleFrequency * 2)
        {
            while (pulseQueue.Count > 5)
                pulseQueue.Dequeue();

            int[] cont = pulseQueue.ToArray();
            float mean = (float)(Mean(cont, cont.Length) + 0.5);
            pulseValue = mean;

            windowCounter = 0;
        }
        else
        {
            windowCounter++;
        }

        // Detect QRS
        int delay = detector.QRSDet(dataByte * 10, 0, sampleCount);
        if (delay != 0)
        {
            int DetectionTime = sampleCount - delay;

            // R-R (in samples)
            long rr = DetectionTime - oldPeak;

            // Calculate R-R (ms) and heart rate instantaneous (bpm)
            Detector.QRS qrs = new Detector.QRS();
            qrs.position = DetectionTime;
            qrs.rr = (short)Detector.GetRR((int)rr, sampleFrequency);
            qrs.bpml = (short)Detector.GetPulse(qrs.rr);

            if (qrs.bpml > 20 && qrs.bpml < 255)
                pulseQueue.Enqueue(qrs.bpml);

            oldPeak = DetectionTime;
        }
    }
    catch (Exception)
    {
        Throw;
    }
}
```



```

/// <summary>
/// Calculate mean of array values.
/// </summary>
/// <param name="value"></param>
/// <param name="nvalues"></param>
/// <returns></returns>
private float Mean(int[] value, int nvalues)
{
    float sum = 0;
    float nval = 0;


    for (int i = 0; i < nvalues; i++)
    {
        sum += value[i];

        if (value[i] > 0)
            nval++;
    }

    if (nval > 0)
        sum /= nval;

    return sum;
}

```

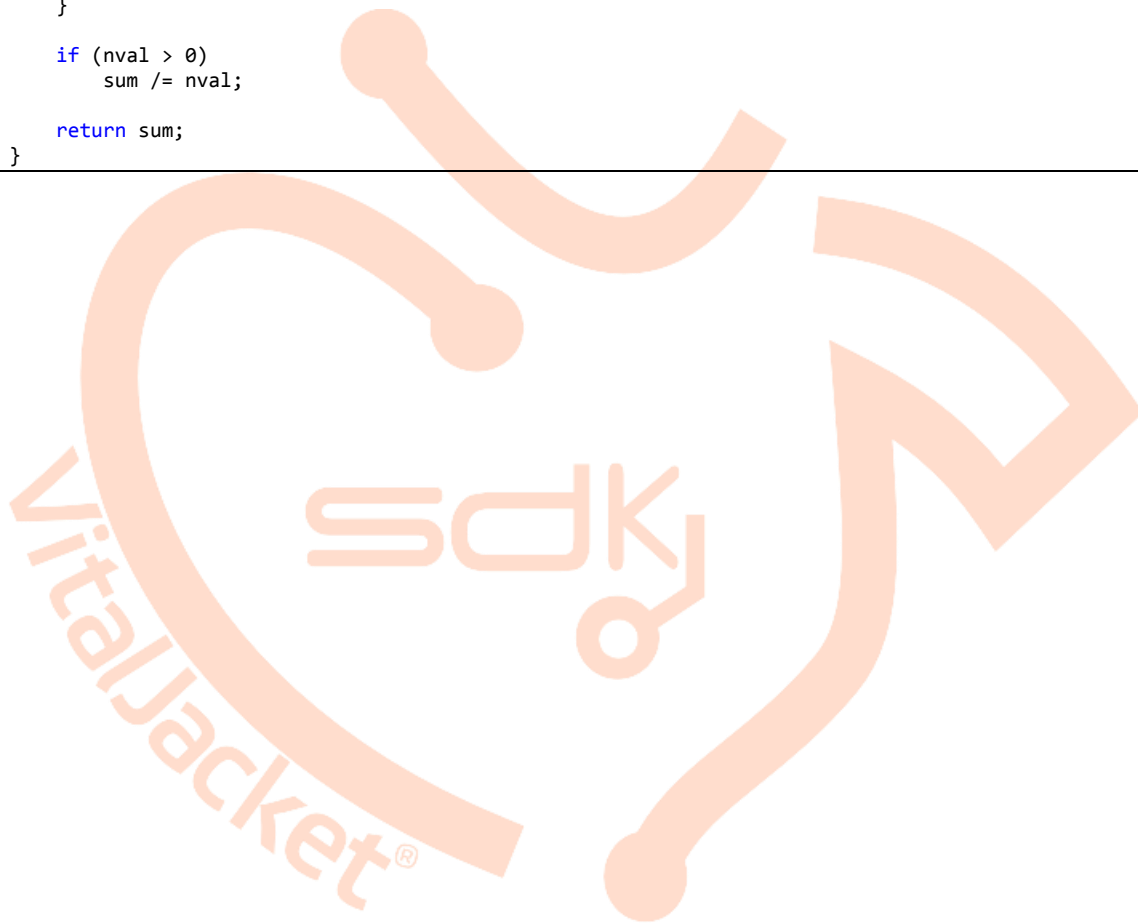

References

- [1] Pan J and Tompkins WJ. A Real-Time QRS Detection Algorithm. IEEE Transactions on Biomedical Engineering 32(3):230-236, 1985
- [2] MIT-BIH Arrhythmia Database: <http://www.physionet.org/physiobank/database/mitdb/>



Control versions

Version	Date	Change log
1.0.02	30-04-2013	<i>Get device Id</i>
		<i>Send radio event to device</i>
1.0.03	19-07-2013	<i>New method to send radio-event to device</i>
		<i>A new too l(InfoExporter.exe) for export data to Excel and Matlab</i>